



# FUJIIRYŌKI

Home use low-frequency wave and potential therapy combination medical device

## FX-14000

ELEDOG 14000

*The best personal care.*



Equipped with "1/f relax waveform" and high performance "prevention from body habituation"

MADE IN JAPAN



## Fluctuation function

Equipped



NEW

[1/f relax] is said to be a rhythm that creates a relaxed state of human body and mind.

Two treatment functions with one unit. Can be used by switching

Equipped with potential therapy and low-frequency wave therapy in a stylish body

## Potential therapy function

In potential therapy a certain amount of potential is given to the body


## Low-frequency wave therapy function

The treatment that treats the affected area by providing a low-frequency current weaker than the skin surface.

## Versatile function, simple operation

### Three waveforms to realize electrification treatment

Equipped with three waveforms for effective treatment at all times. Reduces habituation of the body with electrification treatment with changes.

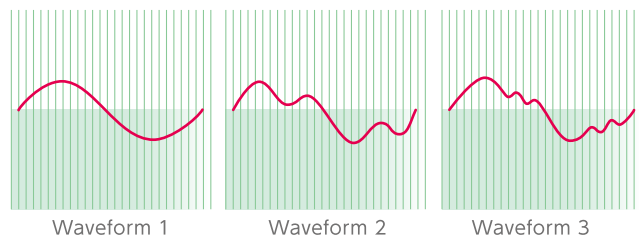

### Equipped with 1/f fluctuation waveform that automatically changes voltage

[1/f relax] is the rhythm that gives relax mood to human body such as the sound of wave, rushing creek. The rhythm of therapy can be changed. By pressing 1/f button. ON/OFF is available with 1/f button.

### It performs optimum and effective setting for user's body condition or physical capacity.

The voltage strength is 4 levels: 1000Vp, 4500Vp, 9000Vp, 14000Vp and 3 waveforms, 4 levels of energizing time is 20 minutes, 40 minutes, 60 minutes, and 8 hours. If you are using for the first time, please start with a low voltage for a short time.

### Detector

Potential therapy device is a treatment with less experience. We use a detector to prevent you from forgetting that you are using the device or noticing the end of treatment.



### 4 courses with various treatment modes (Easy, Fast, Comfortable, Manual)

Course	Easy	Fast	Comfortable	Manual
Target	For beginner	Person who get used to	Prolonged time	For manual operation
Output voltage (Peak value)	Fixed 4500Vp	Fixed 14000Vp	Fixed 1000Vp	Choice 4500Vp/9000Vp 14000Vp/Random
Output wave	Fixed Wave 1	Random (Wave1~3)	Fixed Wave 1	Choice Wave (wave1~3)
Offset rate	Fixed 45:55	Random 45:55/35:65 25:75	Fixed 50:50	Choice 45:55/35:65 25:75
Frequency	Fixed 60Hz	Random 40·60·100Hz	Fixed 60Hz	Choice 40·60·100Hz
Timer	Fixed 20min	Fixed 60min	Fixed 8hours	Choice 20·40·60min

Fast course: (Output waveform, offset ratio, frequency) automatically change by 3 minutes. Same output waveform, offset ratio, frequency may occur.

### New function to prevent habituation.

Equipped with a function for the purpose of prevention from body habituation that is brought by getting same treatment repeatedly. Most effective setting according to each user's physical condition is available. \*Only for manual course.

### Operability

Easy operation of just press the button next to the display while watching LCD.



- **Frequency**  
Select from 40, 60, 100Hz.  
(The number of vibration per second.)
- **Offset ratio**  
Select from 45 : 55, 35 : 65, 25 : 75,  
which prevents habituation by changing the ratio of plus and minus.


Use the product according to your physical condition. Unreasonable use may cause poor physical condition.

**Please note: Touching other people or objects (metallic) while the power is on energized, it may occur a painful stimulation.**

Gentle stimulation for pain or stiffness. Easy to use low-frequency wave function.

### 5 automatic course settings corresponding to the characteristics of the affected area.

5 automatic course of Shoulder/Waist/Knee/Arm & Elbow/EMS+M.Current. Programs according to each affected area are preset.

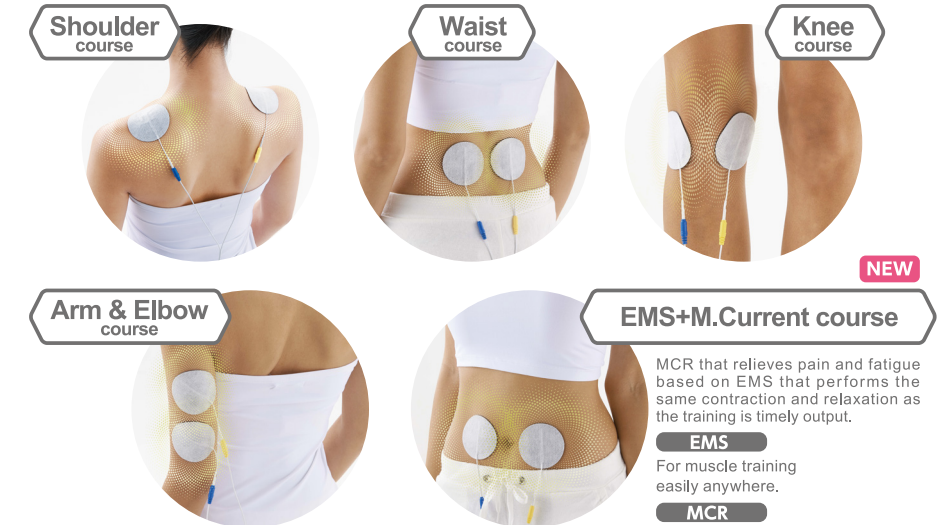

### Various kneading waveform that is output according to set course.

- Wavelet
- Tapping A
- Tapping B
- Knead A
- Knead B
- Press
- Vibration A
- Vibration B

Automatic course is output from 8 kneading waveform of Wavelet/Tapping A, B/Knead A, B/Press/Vibration A, B according to the characteristic of affected area.

※Except for "EMS+M.Current" course

### Potential therapy Purpose of use or effects

### Relief of headache, stiff shoulders, insomnia and chronic constipation. Home use only

• Hypothesis of home use potential therapy device.  
The effect that seems to affect the efficacy or effect is that because it is confirmed that the body hair shakes, stimulates the skin, and that the skin temperature rises when an electric field is applied to humans. It is thought that it stimulates sensory receptors that sense tactile sensation and pressure, and works on blood circulation and body regulation. This "action hypothesis" is a hypothesis of action obtained by an evaluation committee by medical, scientific and engineering experts based on the results of research on literature, etc., not an effect. Source: Japan Home Health Equipment Association "Adequate Advertising and Display Guide for Home Medical Equipment III" 2011 edition

### Low-frequency wave therapy device Purpose of use or effects

Relief of stiff shoulders, prevent paralyzed muscle atrophy and massage effect (Improve blood circulation, recover from fatigue, relief of muscle pain and nerve pain, relief of muscle stiffness, relieve muscle fatigue). Home use only.

Home use low-frequency wave and potential therapy combination medical device

# FX-14000

## ELEDOG 14000

Dimensions(Approximately)/ ■ Main unit: W390 x L340 x H145 (mm)  
 ■ Insulation mat: W1,200 x L2,000 (mm) ■ Energizing sheet: W280 x L595 (mm)  
 ● Weight(Approximately): 7.8 kg ● Power voltage: AC100V(50/60Hz)  
 ● Length of power cord(Approximately): 3m ● Length of the cord of energizing sheet(Approximately): 1.5 m  
 ● Medical device certification number  
 227AKBZX00119000(Managed medical device)  
 ● JAN code: 4951704432122

- |   |   |
|---|---|
| <p><b>■ Potential therapy device</b></p> <ul style="list-style-type: none"> <li>● Power consumption(Approximately): Maximum 20W (When 14,000 Vp)</li> <li>● Output voltage: 1,000 Vp, 4,500 Vp, 9,000 Vp, 14,000 Vp</li> <li>● Timer: 20, 40, 60(min), 8 hours</li> </ul> | <p><b>■ Low-frequency wave therapy device</b></p> <ul style="list-style-type: none"> <li>● Power consumption(Approximately): Maximum 10W</li> <li>● Oscillation frequency(Approximately): Less than 1.2 KHz</li> <li>● Timer: 15 minutes</li> </ul> |
|---|---|



Supplied accessories	① Energizing sheet/1set ② Insulation mat/1set ③ Voltage detector with AA manganese batteries / 1 set ④ Power cord/1unit ⑤ Low-frequency cord/1unit ⑥ Low-frequency pad(4 in 1) / 1 unit
Consumable accessories	① Energizing sheet ⑥ Low-frequency pad
Optional accessories	⑦ Energizing sheet/Long type: W280 x L1,000 (mm)

Easy to operate at a touch of button while watching the screen.

◆ Press "Power" button and display the screen. The screen changes to below screen after 2 seconds.

◆ Background of screen can be changed by pressing "1/f relax" button on the setting screen.

◆ Night mode The screen can be darkened by pressing "Select button" and "OK button" at the same time during treatment.

◆ The sound volume can be decreased by pressing "▲Strong button" and "▼Weak button" at the same time when selecting treatment mode and course.

◆ Select potential therapy or low-frequency therapy with "Select button" and press "OK" button.

◆ Choose favorite course with "select button" and press "OK button".

◆ Press "OK button" after choosing favorite course with "Select button".

◆ Starting the therapy with "▶|| Start/Stop" button.

◆ Starting the therapy with "▶|| Start/Stop" button.

**※ON/OFF the 1/f function with "1/f button".**



Please observe for your safety.

**Danger**

- Do not use it in combination with medical devices such as the following:
- Pacemaker, implantable type defibrillator and other implantable medical electrical devices which are susceptible to electromagnetic interference
- Mountable medical electrical devices such as electrocardiograph.
- Do not use it if you are diagnosed with heart disease and are restricted from daily extreme exercise.

**Warning**

- The following persons must consult a healthcare professional before using this product. ● If you are suffering from spinal fractures, sprains, muscle tears, etc. ● If you are pregnant or have recently given birth ● If you have malignant tumors ● If you have heart diseases ● If you are prescribed rest ● If you have perceptual disorder caused by advanced endovascular disorder due to diabetes etc. ● If you have a body temperature above 38 °C (with a fever) ● Acute inflammatory symptoms (malaise, Chills, blood pressure fluctuations, etc.) ● Those with infectious diseases ● Those who cannot express their intention ● Those who are drunk ● If you are receiving medical treatment, especially when you feel abnormal. ■ The person who can handle the following must be accompanied by a person. ● Those who need nursing care, children ● Those who can not understand the contents of the instruction manual. ■ When any other person uses (friend etc.) accompany them or let them use only after explaining. Please explain to them that they can use it after reading the instruction manual. ■ Do not use for infants and people with physical disabilities. ■ Do not use when drunk or when you have taken sleeping pills. ■ Do not use if the electrode is wet. ■ People around you should not touch you during the treatment. ■ If you have any abnormal feelings in the body stop using immediately.

MEMO

MADE IN JAPAN

Note: specification subject to changes without notice due to product improvements. Some details may differ from actual product.