

## To clean and maintain the product

### 1. Cleaning

Before you clean, please make sure that the switch is in the off location. And the power plug has been pulled out from the socket. Please use dry cloth to clean the product.

### Attention

- Please make sure that do not use benzene or thinner and other aggressive detergents to clean the product. Or the product may fade.

### 2. Maintenance

Make sure that the switch is off, and the power plug has been pulled out from the socket. Keep the product clean and far from the wet environment.

## Technology parameter

Name : Wireless Massage Cushion  
Model : MC162  
Rated Voltage : DC12.6V 2A  
Rated Power : 25.2W  
Battery Capacity : 11.1V 1800mA  
Made in China

Johnson Health Tech. Co., Ltd.  
No. 999, Sec.2, Dongda Rd. Daya Dist., Taichung City, Taiwan,  
R.O.C.



### Correct Disposal of this product


This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

# SYNCA

USER MANUAL

## i-Puffy PLUS

Wireless Massage Cushion  
MC162

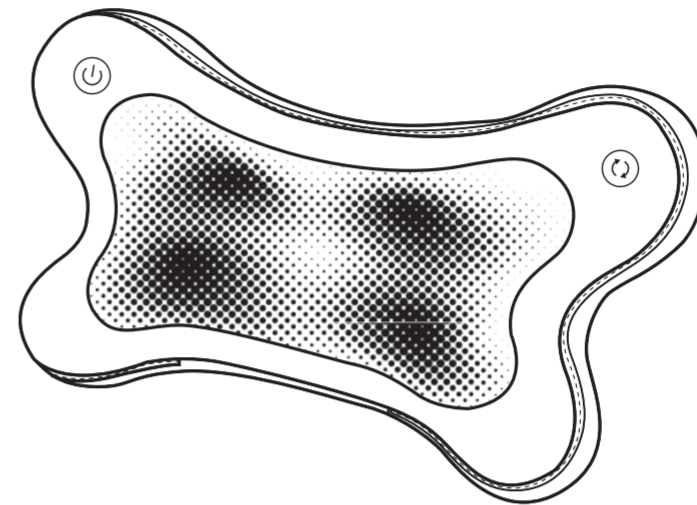

This manual should be kept on hand for easy reference.  
It is your guide to a safe and efficient operation.

## IMPORTANT SAFETY INSTRUCTIONS

When using the product, all instructions and basic precautions should always be followed, including the following:

### DANGER – To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

**WARNING** – For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance. To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not in use and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire. Electric shock, or injury to persons.
- 3) Close supervision is necessary when this product is used by, or near children, invalids, or disabled persons.
- 4) Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this product if it is dropped, has a damaged cord or plug, a malfunction, or it is dropped into water. Return the product to service centre for examination and repair.
- 6) Do not carry this product by its supply cord or use cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair and the like.
- 9) Never drop or insert any object into any openings.
- 10) Do not use outdoors.
- 11) Do not operate where aerosol (Spray) products are being used or where oxygen is administered.
- 12) To disconnect, switch off all controls and remove plug from outlet.
- 13) This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- 14) Children shall not play with the appliance, cleaning and user maintenance shall not be made by children without supervision.
- 15) The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- 16) Do not use massager in close proximity to loose clothing or jewelry.
- 17) Keep long hair away from massager while in use.
- 18) This appliance contains batteries that are only replaceable by skilled persons.
- 19) this appliance contains batteries that are non-replaceable.
- 20) The appliance is only to be used with the power supply unit provided with the appliance.
- 21) The instructions are also available on <https://www.synca-wellness.com>

SAVE THESE INSTRUCTIONS  
This product is intended for household use only

## How to operate

1. Plug in the mains plug and charge a battery.
2. Press "Power Switch" for 3 seconds to start massage and to power off.
3. Press "Function Switch" to change the kneading direction. The preset is clockwise kneading massage(2D massage); Press "Function Switch" to change to anti-clockwise kneading massage(3D massage). Press "Function Switch" and hold for 3 seconds to turn ON / OFF heat function.
4. 20 minutes timer set, which is convenient and safe.

After use : When the product is not in use, it should be unplugged.

The massager is only to be used with the power adapter provided with the appliance.

### Attention

RISK OF FIRE AND ELECTRIC SHOCK.  
FOR DRY LOCATION USE ONLY.  
DO NOT SHORT CIRCUIT.  
DO NOT DISASSEMBLE OR REASSEMBLE THIS DEVICE.  
DO NOT USE IT WITHOUT INSTRUCTION.

## Massage position



## Name and function of the machine parts

